

# How Does HotDog Compare?

|                                        | HotDog            | Hot Air           | Hot Water         |
|----------------------------------------|-------------------|-------------------|-------------------|
| Effective Warming                      | ✓                 | ✓                 |                   |
| Even Heat Distribution                 | ✓                 |                   |                   |
| Low Ongoing Costs                      | ✓                 |                   |                   |
| Durable Waterproof Blankets            | ✓                 |                   |                   |
| Easy to Clean                          | ✓                 |                   |                   |
| Won't Dry EKG Lead or Mucous Membranes | ✓                 |                   | ✓                 |
| Maintenance-free Controller            | ✓                 |                   |                   |
| Quiet Operation                        | ✓                 |                   | ✓                 |
| Fast Warm-up Time                      | ✓                 | ✓                 |                   |
| Small and Light for Ease of Use        | ✓                 |                   |                   |
| No Blowing Air on Surgical Site        | ✓                 |                   | ✓                 |
| No Hot Air Blowing on You              | ✓                 |                   | ✓                 |
| Cable Easily Fits through Cage Slats   | ✓                 |                   | ✓                 |
| Environmentally Friendly               | ✓                 |                   | ✓                 |
| Energy Efficient                       | Average 100 watts | Average 900 watts | Average 800 watts |



We Sell Direct.

Over 30 years ago, I invented Bair Hugger® warming. Now, to address your concerns about air and water systems, I have developed the first conductive fabric warming for veterinary use.

HotDog Patient Warming is great for cats, dogs, exotics and more. I am confident it will help you provide the most effective, safe, and economical solution to meet all of your patient’s warming needs.



Warm regards,  
Scott Augustine, M.D.



“The HotDog Patient Warming System provides superior means for perioperative thermoregulation.”

-Iowa State University clinical study



Contact Us Today

(866) 484-3505

cs@aug surg.com

www.vetwarming.com



Follow us on Social Media

@vetwarming



©2025 Augustine Surgical, Inc. All Rights Reserved.  
Bair Hugger is a registered trademark of Solventum.

6581 City West Parkway,  
Eden Prairie, MN 55344



Featuring  
ThermAssure™  
Conductive Fabric





# No Hot Air... No Water... Just Safe, Even Warmth

HotDog Patient Warming represents a significant breakthrough in technology and features our exclusive ThermAssure fabric.

ThermAssure's conductive polymer responds to low-voltage, direct current by becoming uniformly warm without hot spots.

The System includes many built-in safety features.



Clinically proven to be more effective than forced-air and circulating water systems.

### Controllers:

- Two Controller options: run one or up to four Blankets
- Touchscreen user interface for easy operation
- Lightweight
- Maintenance-free
- Quiet operation
- Free-standing or IV pole mountable (VESA)

### Blankets:

- Conductive fabric technology
- No wires as a heat source
- X-ray, translucent fabric
- Use over, under, and around patients
- Great for dentals, surgery, and recovery
- Durable, reusable, and cleanable

## HotDog Veterinary Warming System



| Controller | Description            | Primary Use                     |
|------------|------------------------|---------------------------------|
| WC71V      | Single Port Controller | Runs 1 Blanket                  |
| WC77V      | Multiport Controller   | Runs up to 4 Blankets           |
| Blanket    | Description            | Primary Patient Use             |
| V101       | X-Small 10" x 12"      | Mice and Rats in Research       |
| V102       | Small 10" x 23"        | Exotic or Pocket Pets           |
| V103       | Medium 16" x 22"       | Average Sized Cat or Small Dogs |
| V104       | Large 22" x 31"        | Average Sized Dogs              |
| V106       | X-Large 27" x 47"      | Large Dogs                      |
| Accessory  | Description            | Primary Use                     |
| A101       | Cable 13' (standard)   | Connects Blanket to Controller  |
| A101-6000  | Cable 20'              | Connects Blanket to Controller  |
| V501       | Exam Table Cover       | Used on most surfaces           |
| V600       | Pressure-Reduction Pad | Reduces risk of pressure injury |


## Pressure-Reduction Pad



- Reduce the risk of pressure injury
- Use over a hard surface, surgery table, cage, or dental grate
- Easy to clean
- 20" x 31"



## Reduce Pressure in Risk Zones

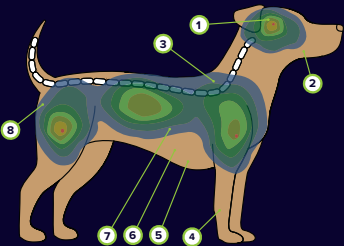

1. Cranium
2. Mandible
3. Scapula
4. Carpal Bones
5. Sternum
6. Cartilaginous Part of Rib
7. Ribs
8. Ischium

## Less Stress™ Exam Table Covers

- Stop paws from slipping
- Pets feel safer
- Cover cold table surfaces
- Easy-to-clean Polyethylene material



Adhesive backing leaves no residue when removed.

Exam Table Covers can be used on most surfaces, including formica and stainless steel. The convenient 20" x 40" size fits most exam tables. It's great for scales, too!